

mindray

Accessories for operating tables

www.mindray.com

P/N:ENG-Accessories for operating tables-210210x52P-20200110
©2020 Shenzhen Mindray Bio-Medical Electronics Co.,Ltd. All rights reserved. Specifications subject to changes without prior notice.

mindray
healthcare within reach



Accessories for Operating Tables

Accessories for operating tables from Mindray are characterised by these functionality, flexibility and high reliability. Developed in close collaboration with end users, the accessories fulfill all kinds of surgery positioning demands. They place the patients in an optimum way that offer benefits both for surgeons and patients.

We use high quality stainless steel materials in our operating tables and accessories. In order to get the surgeons access to C arm, all the pads or plates are radiolucent for X-ray check.

To avoid discomfort in sensitive areas of the patients, we provide gel material pads for patient positioning which are similar as human tissue.



Supine position



Prone position



Lateral position



Lithotomy position



mindray
healthcare within reach



01-12 General Accessories

01



13-18 Head Accessories

13



19-22 Leg Accessories

19



23-24 Back Accessories

23



25-28 Accessories for Gynecology & Urology

25



29-36 Orthopedic Accessories

29



37-38 Accessories For Controlling
And Easier Application

37




39-42 Comfort solution Accessories

39

HyBase 8300/8500

Electro-Hydraulic Operating Table

Engineered for increasingly complex surgical scenes, HyBase 8300/8500 comes with extreme weight capacity up to 460 kg guaranteeing highest safety and stability in any position. With longitudinal slide up to 480 mm, it is ideally suitable for cardiovascular surgery with C-arm. The modular design tabletop enables it to be tailored for various surgical discipline needs.

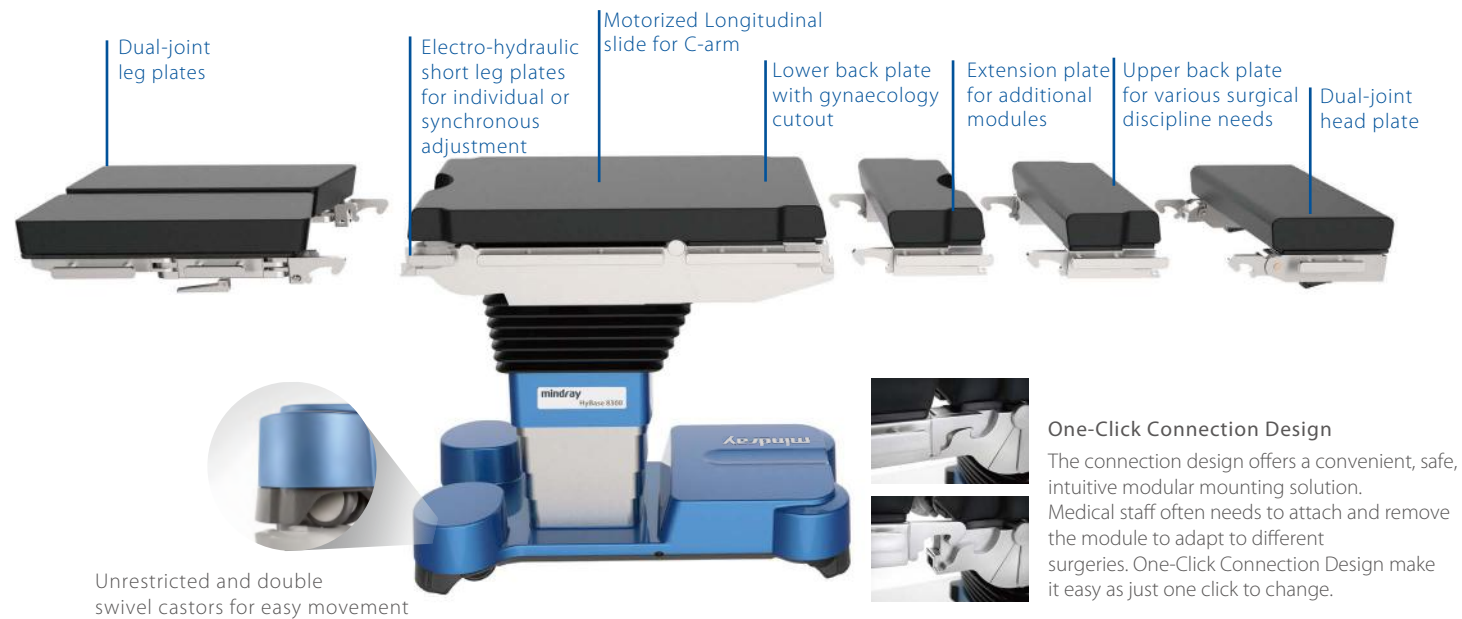

- Longitudinal slide up to 480 mm for free access to C-arm
- Extreme weight load capacity of 460 kg with highest safety and stability
- Modular design tabletop tailored for various surgical discipline needs
- Intelligent anti-collision system to avoid table components collision during movement
- Multi-layer decompression pad with waterproof and seamless design
- 575 mm lowest table position for Neurosurgery
- IR wireless remote control system for free movements in OR

HyBase 6100

Electro-Hydraulic Operating Table

Providing a better clinical use's scenario, HyBase 6100 is designed with a much lower table height of 498mm specially for neurosurgery and minimally invasive surgery, a more user friendly base cover, as well as a more convenient wireless remote control system. Moreover, with longitudinal slide up to 320 mm, the table is ideally suitable for cardiovascular surgery or spinal surgery with C-arm. All in all, HyBase 6100 maximizes the flexibility of operating table's versatility.



- Longitudinal slide up to 320 mm for free access to C-arm
- Max. lifting capacity up to 380kg with highest safety and stability
- Max. safe working load 275kg in normal position
- Multi-layer decompression pad with waterproof and seamless design
- 498 mm lowest table position for Neurosurgery and Minimally Invasive Surgery
- Wireless remote control system
- Foot pump Manual Override System


HyBase 3000

Electro-Hydraulic Operating Table

Proper annotation of "cost-effective". with most acceptable price, gaining the most functional Electro-hydraulic table. No matter what principle surgery was, it can always find a way to meet your requirement.



Foot break provides Maximum safety and easy use.

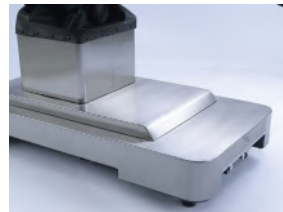

Battery power can support 50-80 operations.
Backup control keys in case of hand control failure.
Optional stainless steel cover.

- Multi-layer decompression pad with waterproof and seamless design
- Longitudinal slide up to 300 mm for free access to C-arm
- Max. lifting capacity up to 230kg
- Safe working load up to 185kg in normal position
- Battery power for about one week time
- Mechanical brake system for optimum stability
- Return to original position with one button
- Flex & Reflex position with one button
- Radiolucent table plate for free C-arm access
- Three control panels including remote panel, backup control & foot switch

UniBase 30

Electro-Hydraulic Operating Table

During the surgery, a good operating table is the cornerstone of success. Engineered for increasingly complex surgical scenes, UniBase 30 needs least effort and time to position the patient in a safe, comfortable and efficient way, providing an economic solution without straining on your hospital's budget.



Versatile & Ergonomic Design



Detachable Head Plate



Adjustable Body Elevator(optional)



Pedal Floor Lock

- Weight capacity of 185 Kg in normal position
- UniBase30, simple good operating table:
- Longitudinal shift (optional) with 300mm
- Adjustment body elevator (optional) with 120mm
- Battery power for about one week time
- Mechanical brake system for optimum stability
- Detachable & split leg plate for lithotomy position

General Accessories

The Mindray table accessory portfolio provides optimal patient positioning and ideal working conditions for the medical staff. With robust design, ease-to-use and ergonomic functionality, Mindray offers standard table accessories that meets the requirements of nearly every surgical discipline.



General Accessories



0092-30-102657

Light quadrate clamp

Accessory Name:
Universal clamp (Light, one-piece, quadrate)

Application:
To attach to the side rail for connecting accessories

Compatible with:



0092-30-102655

Heavy quadrate clamp

Accessory Name:
Universal clamp (Heavy, one-piece, quadrate)

Application:
To attach to the side rail for connecting accessories with the table

Compatible with:



0092-30-102658

Light radial clamp


Accessory Name:
Universal clamp (Light, one-piece, radial)

Application:
To attach to the side rail for connecting accessories with the table

Compatible with:



General Accessories



Heavy radial clamp

Accessory Name:
Universal clamp (Heavy, one-piece, radial)

Application:
To attach to the side rail for connecting accessories with the table

0092-30-102656


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



Strap clamps, a pair

Accessory Name:
Strap clamps, a pair

Application:
Can fix the body strap at a specified location on the side rail

115-038243-00


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



Infusion holder with clamp

Accessory Name:
Infusion pole with one piece of heavy radial clamp

Application:
To hang the infusion bag in the procedure

0092-30-102643

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

General Accessories



Light Ana. frame with clamp

Accessory Name:
Anesthesia Frame, one-piece with light quadrate clamp

Application:
To hang the anesthesia screen during the surgery

0092-30-102645


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



U-Shaped Ana. frame w/ clamps

Accessory Name:
"U" shape Anesthesia Frame with clamps

Application:
To hang the anesthesia screen during the surgery


115-005689-00

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase



Heavy Ana. frame with clamp

Accessory Name:
Anesthesia Frame , "L" shape design with two adjustable tubes, one-piece with clamp

Application:
To hang the anesthesia screen during the surgery

0092-30-102654

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

General Accessories



Light quadrate clamp

Accessory Name:
Arm board with a rotation function in a horizontal position. One piece, with an integrated clamp, a pad and two pieces of fastening belt

Application:
To support and fasten the arm of the patient to avoid unexpected moving

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

0092-30-102646



Heavy arm board (PU/MFP) with clamp

Accessory Name:
Arm board with an up and down and a rotation function. One piece, with one heavy radial universal clamp, a pad and two pieces of fastening belt (SFC/PU)

Application:
To support and fasten the arm of the patient to avoid unexpected moving

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-016721-00 0092-30-102634

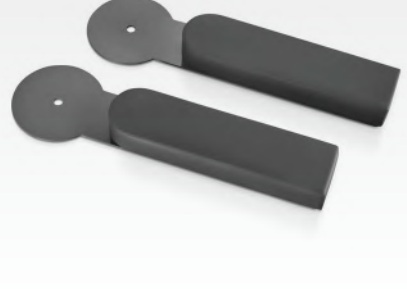

Carbon fiber arm board, a pair

Accessory Name:
Carbon fibre arm board. A pair.


Application:
A carbon fibre arm board which includes a memorised foam pad. This can work with the 1180 cardiac operating plate which provides support for the patient's upper extremities in a supine position

Compatible with:

8500/8300
HyBase

115-030153-00

General Accessories



Double arm board with clamp

Accessory Name:
Double layer arm support with one piece of heavy radial clamp for lateral position, with the pad

Application:
To support and fasten both of the arms of the patient as a lateral position to avoid unexpected movement

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-005693-00



Arm board for lateral position w/ clamp

Accessory Name:
Lateral position arm support with onepiece of radial universal clamp, with the pad, one-piece

Application:
To support and fasten the upper arm of the patient in lateral position

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

0092-30-102648



Arm posturing device

Accessory Name:
Arm posturing device

Application:
During surgeries, it is used to allow support for patient's arm to effectively protect patient's brachial nerves

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

115-038996-00

General Accessories



115-007051-00

Suspended arm strap

Accessory Name:
Suspended arm strap

Application:
To attach to anesthesia frame for hanging the hand of the patient


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-007053-00

Wristlet with clamp

Accessory Name:
Wristlet with one-piece of quadrate universal clamp

Application:
To attach to the side rail for fastening the hand of the patient

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



0092-30-102640

Body strap

Accessory Name:
Body strap, one-piece

Application:
To fasten the body of the patient to the table during the surgery

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

General Accessories



0092-30-102669 0092-30-102638

Heavy body support with clamp, pubis/sacrum

Accessory Name:
Heavy body support with clamp, pubis
Heavy body support with clamp, sacrum

Application:
To support the body of the patient for avoiding unexpected movement

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



0092-30-102650

Light body support w/ clamps, a pair

Accessory Name:
Two pieces of body support with a pair of quadrate universal clamps and the pad

Application:
To support the body of the patient for avoiding unexpected movement

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



0092-30-102637 0092-30-102671

Heavy shoulder support with clamp, left/right

Accessory Name:
Left/right shoulder support with pad (one piece, Heavy version with one piece of heavy quadrateuniversal clamp and pad)

Application:
To support the shoulder of the patient for avoiding unexpected movement

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

General Accessories



Light shoulder support w/ clamps, a pair
Accessory Name:
Shoulder support (left and right, pair with quadrate universal clamps and pad)

Application:
To support the shoulder of the patient for avoiding unexpected movement

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

0092-30-102649



Cassette fixer with clamp
Accessory Name:
X-ray cassette fixer (compatible with multi-size)

Application:
To fix the X-ray cassette under the table

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

0092-30-102641



Pad for intervertebral disc operations
Accessory Name:
Intervertebral disc operation pad (Memorise foam, one piece)

Application:
To support the body in intervertebral operating position

Compatible with:

8500/8300
HyBase


6100
HyBase

3000
HyBase

30
UniBase

0092-30-102639

General Accessories



Bow frame for spine operation
Foam head block stand for bow frame(optional)
Pad set for bow frame,a pair
Accessory Name:
Bow frame for spine operation,Foam head block (optional),
Pad set for bow frame,a pair(optional)

Application:
Bow frame is adjustable for length and width according to different patients.
Foam head block support the head in prone position
Pad set give a comfortable support for the chest.

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-052219-00 / 115-053656-00 / 115-053654-00



Bow frame cart
Accessory Name:
Bow frame cart, four universal castors

Application:
It is convenient for the bow frame to store and move

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-052221-00

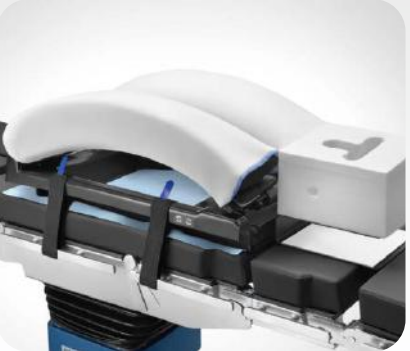

Disposable covers for bow frame(6 pcs)
Accessory Name:
covers for bow frame,disposable

Application:
To prevent the pad getting dirty during spine operation

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-052223-00

General Accessories



115-033429-00

Side rail extension

Accessory Name:

Side rail extension

Application:

Can extend the side rail toward the two ends of the operating table, to provide more space for installation of other accessories, such as leg support and foot plate

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-039082-00

Side rail extension

Accessory Name:

Side rail extension (HySRE-2)

Application:

Can extend the side rail towards the two sides of the operating table, to provide more space for installation of other accessories

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-030814-00

Table width extender with clamps, a pair

Accessory Name:

Table width extender with clamps, a pair

Application:

During surgeries, it is used to widen the operating table top and provide comfortable lying for larger-sized patients

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



Head, Back and Leg Accessories

Modular design table top is tailored for various surgical discipline needs, which also can low down the cost of the hospital as they only need to buy the specific module for more surgeries not one more table.



Head Accessories



Head plate
Head plate, CF
Accessory Name:

Head plate, single-joint, carbon fiber (optional), with special foam core pad.

Application:

During surgeries, it is used to provide support to the patient's head

Compatible with:



115-028574-00 115-028576-00



2-Joint head plate
2-Joint head plate, CF

Accessory Name:

Head plate, dual-joint, carbon fiber(optional), with special foam core pad

Application:

During surgeries, it is used to provide support to the patient's head with a double-joint design

Compatible with:



115-028703-00 115-028705-00



Back plate for shoulder operation
with helmet

Accessory Name:

Demountable shoulder support module with helmet to support patient's head.

Application:

Provide a stable support for the shoulder operation by beach-chair position

Compatible with:



115-052229-00

Head Accessories



115-012835-00

Head rest, Type 2

Accessory Name:
Adjustable horseshoe-shaped head rest with two installation bars

Application:
To support the head of the patient providing a stable support for Neurosurgery operation.You should select the adaptor for head position simultaneously


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-059309-00

Head rest

Accessory Name:
Head rest, new memorized foam pad

Application:
A solution for supporting the head safely and comfortably.You should select the adaptor for head position simultaneously


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



0092-30-102642

Adaptor for head position

Accessory Name:
Neurosurgical head support adapter

Application:
To connect the adjustable horseshoe, headrest,headrest system.


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

Head Accessories



115-034319-00

Universal adaptor for head position

Accessory Name:
Universal adaptor for head position

Application:
To connect the adjustable horseshoe, headrest,headrest system.

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-007045-00

Ophthalmic operating plate

Accessory Name:
Ophthalmic operating head plate (compatible with the adaptor of head position)

Application:
To support the head of the patient and supprt the arm of the surgeon during the surgery


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-007045-00

Wrist Support

Accessory Name:
Wrist Support ,foam pad, used with head rest

Application:
Supporting the wrist stably to reduce the surgeon's fatigue

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

Head Accessories




Guiding roller for cervix traction

Accessory Name:
Guiding roller for cervix traction,used with head rest

Application:
It is used for the cervical vertebra traction o of the spine operation

Compatible with:
8500/8300 HyBase 6100 HyBase 3000 HyBase 30 UniBase

DJ15-0101-101




Aluminum headrest system

Accessory Name:
Aluminum headrest system (Including base unit, swivel adaptor and skull clamp; not including skull pins)

Application:
For supine, prone, lateral positions, please select the adaptor for head position(PN:0092-30-102642 or 115-034319-00); For sitting positions, please select crossbar adaptor (PN: DJ15-0101-101); only allow to sale with the table together as a package

Compatible with:
8500/8300 HyBase 6100 HyBase 3000 HyBase 30 UniBase

DJ15-0101-001




Radiolucent headrest system

Accessory Name:
Radiolucent headrest system (Including base unit, swivel adaptor, transitional member and skull clamp; not including skull pins)

Application:
Radiolucent for C-arm access, only allow to sale with the table together as a package

Compatible with:
8500/8300 HyBase 6100 HyBase 3000 HyBase 30 UniBase

DJ15-0101-002



Adult horseshoe headrest w.Ext Pedia horseshoe headrest w.Ext

Accessory Name:
Adult horseshoe headrest w.Ext
Pedia horseshoe headrest w.Ext

Application:
A solution for supporting the head more safely and comfortably,used with headrest system

Compatible with:
8500/8300 HyBase 6100 HyBase 3000 HyBase 30 UniBase

DJ15-0101-301 DJ15-0101-301




Crossbar adaptor with clamps

Accessory Name:
Crossbar adaptor with clamps for headrest system in sitting positions

Application:
Only allow to sale with the table together as a package;

Compatible with:
8500/8300 HyBase 6100 HyBase 3000 HyBase 30 UniBase

DJ15-0101-101



Titanium skull pin, pediatric (3pcs); SS skull pins, pediatric (3 pcs) Titanium skull pin, adult (3pcs); SS skull pins, adult(3 pcs)


Accessory Name:
Titanium skull pin, pediatric (3pcs);
Stainless steel skull pins, pediatric(3 pcs)
Titanium skull pin, adult(3pcs);
Stainless steel skull pins, adult(3 pcs)

Application:
Only allow to sale with the table together as a package;

Compatible with:
8500/8300 HyBase 6100 HyBase 3000 HyBase 30 UniBase

DJ15-0101-201 DJ15-0101-203
DJ15-0101-202 DJ15-0101-204

Leg Accessories



**1-joint leg plate, a pair ;
1-joint leg plate, a pair, CF**


Accessory Name:
Leg plate, single-joint, carbon fiber (optional), with special foam core pad, a pair

Application:
During surgeries, it is used to provide support to the patient's leg

Compatible with:

8500/8300
HyBase

115-028580-00 115-028578-00



**2-joint leg plate, a pair ;
2-joint leg plate, a pair, CF**


Accessory Name:
Leg plate, dual-joint abduction, carbon fiber (optional), with special foam core pad, a pair

Application:
To support the legs of the patient during the surgery. During surgeries, it is used to provide support to the patient's leg. A double-joint design, helps doctors get closer to the surgical site.

Compatible with:

8500/8300
HyBase

115-028582-00 115-028584-00



**4-parts leg plate;
4-parts leg plate, CF**

Accessory Name:
Leg plate, 4 sections, with special foam core pad carbon fiber (optional), a pair

Application:
Can achieve the abductive function of normal leg plates and horizontal abductive function of 2-joint leg plates; the second joint leg plate can be folded upward and downward, abducted and can extend toward the leg end, and is detachable; include a mattress made of memory foam. It is able to assist in achieving kneeling position, lithotomy position and supine position.

Compatible with:

8500/8300
HyBase

115-033426-00 115-038559-00

Leg Accessories



**One-section leg plate, hook ver;
One-section leg plate, hook, CF**

Accessory Name:
One-section leg plate, hook ver
One-section leg plate, hook, CF

Application:
During surgeries, it is an one-section plate used to support the patient's legs.

115-059309-00

Compatible with:

8500/8300
HyBase

115-035792-00 115-037739-00



**One-section leg plate,
PU/memory foam pad**

Accessory Name:
One section leg plate for PU or SFC (Special Form Core) mattress


Application:
To support the leg with one-section plate during the surgery

Compatible with:

6100
HyBase

3000
HyBase

115-005683-00 115-005685-00



Transfer board, hook version

Accessory Name:
Transfer board, hook version lithotomy position.


Application:
During surgeries, it is used to temporarily support the patient's lower extremities when positioning the patient or transferring the patient, and to optimize the pre-operation and post-operation workflows at lithotomy position.

Compatible with:

8500/8300
HyBase

115-033425-00

Leg Accessories



115-012841-00

Thigh holder with clamps, a pair

Accessory Name:
Fixation for femur, a pair, with two radial universal clamps

Application:
To support the femur of the patient during the surgery


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-012843-00

Shank holder with clamps, a pair

Accessory Name:
Fixation for lower leg, a pair, with two radial universal clamps

Application:
To support the lower leg of the patient during the surgery


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



0092-30-102636 0092-30-102667

Foot plate with clamp, left/right

Accessory Name:
Left/right foot support with one piece of heavy quadrate universal clamp and pad

Application:
To support the foot of the patient for avoiding unexpected movement

Compatible with:


8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

Leg Accessories



115-033430-00

Leg restraint cuff, thigh, a pair

Accessory Name:
Leg restraint cuff, thigh, a pair

Application:
Used to secure the patient's thigh.

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



115-033431-00

Leg restraint cuff, shank, a pair

Accessory Name:
Leg restraint cuff, shank, a pair

Application:
Used to secure the patient's shank.

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

Back Accessories



115-033428-00 115-038564-00


**Short upper back plate ;
Short upper back plate, CF**

Accessory Name:
Short upper back plate
Short upper back plate , CF

Application:
During Surgeries, provide support to the patient's back

Compatible with:

8500/8300
HyBase



115-029430-00 115-029432-00


Extension plate,w/ MFP,CF

Accessory Name:
Extension plate, carbon fiber (optional), with gynaecological cutout, with special foam core pad, length of 250mm.

Application:
With gynecological cutout, can extend the head end or the foot end of the operating table top up to 250mm, to provide more perspective space or more standing space for doctors. Applicable to trauma, spine, gynecology, or urology. Includes a optional memorised foam pad.

Compatible with:

8500/8300
HyBase



115-028586-00 115-028588-00

**Upper back plate;
Upper back plate, CF**


Accessory Name:
Upper back plate, carbon fiber (optional),, with special foam core pad.

Application:
During surgeries, provides support to the patient's back.

Compatible with:

8500/8300
HyBase

Back Accessories



115-029434-00

Cardiac operating plate,CF

Accessory Name:
Cardiac operating plate, carbon fiber, with special foam core pad, length of 1,180mm.

Application:
Cardiac operating plate, can provide 360° X-ray radiolucency, extend the foot end of the table top for 1180mm, allows for cardiovascular surgeries and spine surgeries. Includes a memorised foam pad and stand bars.

Compatible with:

8500/8300
HyBase



115-030151-00


Accessory adaptor,CF

Accessory Name:
Accessory adaptor for Cardiac Operating Plate, carbon fibre, with side rails.

Application:
Provide an additional side rail for the 1180 cardiac operating plate (EU Standard), can be used to install accessories such as arm boards and shoulder supports.

Compatible with:

8500/8300
HyBase

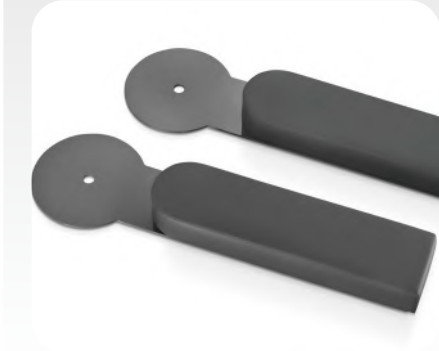

115-030153-00

Carbon fiber arm board, a pair

Accessory Name:
Carbon fiber arm board. A pair.

Application:
A carbon fiber arm board which includes a memorised foam pad. This can work with the 1180 cardiac operating plate which provides support for the patient's upper extremities in a supine position

Compatible with:

8500/8300
HyBase

Accessories for Gynecology & Urology

For Gynecology and Urology Surgery, Mindray provide a better workflow design and ergonomic working space for the patient positioning. The lightweight transfer board will help a lot after positioning the leg with less efforts to remove it. The leg plate folded design helps avoiding the disassembly of the leg plate. And the u shape design on base cover forms a larger space for surgeons while doing the operations. The cutout of the backplate also enable the patient reverse position for Gynecology and Urology Surgery.



Accessories for Gynecology & Urology



115-062765-00

Stirrup (HyStirrup-3), a pair

Accessory Name:
Leg support (Stirrup type), a pair, with clamp

Application:
To support and fasten the leg of the patient

Compatible with:



0092-30-102635

Heavy leg support with clamp

Accessory Name:
One-piece leg support (heavy version). Includes one heavy radial universal clamp, a pad and one piece of fastening belt

Application:
To support and fasten the patient's leg

Compatible with:



0092-30-102647

Light leg support with clamps, a pair

Accessory Name:
Light leg support. A pair with radial universal clamps

Application:
To support and fasten the patient's leg

Compatible with:



Accessories for Gynecology & Urology

A pair of black stirrups with straps, designed for supporting a patient's legs during surgery.

Stirrup (HyStirrup-2)

Accessory Name:
Leg support (Stirrup type), a pair, with clamp without lift assisted

Application:
To support and fasten the patient's leg

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-062768-00

A pair of blue leg holders with loops and clamps, designed for supporting a patient's legs during surgery.

Leg holder with loops, with clamps, a pair

Accessory Name:
Leg Holders. A pair with loops and radial universal clamps.

Application:
To support and fasten the patient's leg.

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-007043-00

A black rectal positioning device with a clamp, designed for supporting a patient's legs during surgery.

Rectal positioning device

Accessory Name:
Rectal positioning device

Application:
To support the coxa and knees of the patient in rectal surgery

Compatible with:

6100
HyBase

3000
HyBase

115-007049-00

Accessories for Gynecology & Urology

A stainless steel drainage bowl with a clamp, designed for collecting liquid during surgery.

Drainage bowl for Gynecology & Urology

Accessory Name:
Drainage bowl for Gynecology & Urology

Application:
To hang the infusion bag in the procedure

Compatible with:

3000
HyBase

115-013018-00

A stainless steel drainage bowl with a clamp, designed for collecting liquid during surgery.

Drainage bowl with clamp, universal

Accessory Name:
Universal drainage bowl

Application:
To drain the liquid in lithotomy position

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-034318-00

A stainless steel drainage bowl with a clamp, designed for collecting liquid during surgery.

Drainage bowl, HyBase 6100


Accessory Name:
Drainage bowl for Gynecology & Urology

Application:
To drain the liquid in lithotomy position

Compatible with:

6100
HyBase

115-005697-00

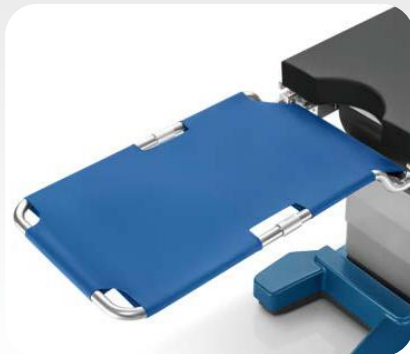

115-052205-00

Transfer board

Accessory Name:

Transfer board,light leg plate

Application:

Support the legs temporarily during lithotomy position

Compatible with:

6100
HyBase



115-033425-00

Transfer board, hook version

Accessory Name:

Transfer board, hook version lithotomy position.

Application:

During surgeries, it is used to temporarily support the patient's lower extremities when positioning the patient or transferring the patient, and to optimize the pre-operation and post-operation workflows at lithotomy position.

Compatible with:

8500/8300
HyBase



115-047824-00

Lateral positioner

Accessory Name:

Lateral positioner,foam pad

Application:

Support the legs during lateral position

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase



Orthopedic Accessories



Orthopedics Extension



Tibia Extension Assembly



Lateral Femoral Extension Assembly



115-066779-00 115-066779-00

Ortho extension w/ cart,CF,basic
Ortho extension w/ cart,CF,comfort

Accessory Name:

Ortho extension device, carbon ber, including the cart

Application:

During the orthopedic traction surgery, it is used to pull the femur to meet the need of x-rays; the device is made of highly-transmissible carbon ber materials to improve image qualities.

Compatible with:

8500/8300
HyBase



115-066781-00 115-066781-00

Ortho extension w/ cart,IV,basic
Ortho extension w/ cart,IV,comfort

Accessory Name:

Ortho extension device, hook version, stainless steel, including the cart

Application:

To provide leg traction during orthopedic surgery

Compatible with:

8500/8300
HyBase



115-066777-00 115-066777-00

Ortho extension w/ cart,III,basic
Ortho extension w/ cart,III,comfort

Accessory Name:

Orthopedic extension device, including the cart

Application:

To provide leg traction during orthopedic surgery

Compatible with:

6100 3000
HyBase

Orthopedic Accessories



115-033410-00 115-066775-00

Ortho extension w/ cart,Uni,basic
Ortho extension w/ cart,Uni,comfort

Accessory Name:
Ortho extension w/ cart, UniBase

Application:
To provide leg traction during orthopedic surgery

Compatible with:



115-007055-00

Tibia & fibula extension assembly

Accessory Name:
Tibia extension assembly

Application:
To provide tibial and fibula traction

Compatible with:




115-007057-00

Lateral femur extension assembly

Accessory Name:
Lateral femoral extension assembly, used with Ortho extension device

Application:
To provide femur traction in lateral position

Compatible with:

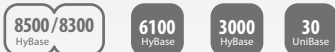


115-047827-00 / 115-047828-00

Leg holder for Ortho extension device,SS
Leg holder for Ortho extension device,CF

Accessory Name:
Support the leg,CF (optional),used with Ortho extension device

Application:
Support the injured leg during the orthopedic operation

Compatible with:

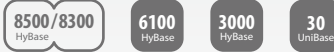

115-011817-00

Pelvis width extender, a pair

Accessory Name:
Triangle support board, 2 pieces.

Application:
It is used to widen the pelvis plate.

Compatible with:




115-033422-00

Upper arm extension device with clamps

Accessory Name:
Upper arm extension device

Application:
During surgeries, it is used to fix and pull patient's forearms; consists of auxillary pillar, elbow joint pillar, palm & finger fixer and tractor

Compatible with:

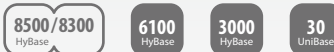

Orthopedic Accessories




115-033421-00

Forearm extension device with clamps

Accessory Name:
Forearm arm extension device

Application:
During surgeries, it is used to fix and pull patient's forearms; consists of auxillary pillar, elbow joint pillar, palm & finger fixer and tractor

Compatible with:

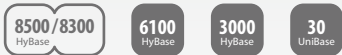


115-033415-00

Upper arm operating table

Accessory Name:
Upper arm operating table

Application:
To support the upper arm of the patient.

Compatible with:

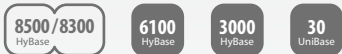

115-007047-00

Arm operating table

Accessory Name:
Arm operating table

Application:
To support the arm of the patient and the hand of the surgeon

Compatible with:



Orthopedic Accessories



115-052229-00 / 115-059315-00

Back plate for shoulder operation with helmet Body strap for shoulder operation

Accessory Name:
Demountable shoulder support module with helmet to support patient's head
Body strap(optional) used with back plate for shoulder operation

Application:
Provide a stable support for the shoulder operation by beach-chair position

Compatible with:



115-033413-00

Arthroscopic leg holder with clamp

Accessory Name:
Arthroscopic leg holder with clamp

Application:
To fasten and support the knee of the patient

Compatible with:



115-033416-00

Meniscus positioning device with clamp


Accessory Name:
Meniscus positioning device with clamp

Application:
To provide support for patient's thigh or foot during knee joint surgeries

Compatible with:



Accessories for Controlling and Easier Application



**IR hand control, LCD, slide;
IR hand control, LCD**

Accessory Name:
IR hand control, with LCD, with/without longitudinal slide function, with wireless charger.

Application:
Can adjust and control movements of the operating table through the infrared remote control system; with LCD. Wireless charging and easy to clean.

Compatible with:

8500/8300
HyBase

115-028600-00 115-028598-00



**Corded hand control, LCD, slide;
Corded hand control, LCD**


Accessory Name:
Corded hand control, with LCD, with/without longitudinal slide function

Application:
Adjust and control movements of the operating table; with LCD

Compatible with:

8500/8300
HyBase

115-028596-00 115-028594-00



Foot switch

Accessory Name:
Foot switch

Application:
To control the table movement


Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

0092-30-102644



Accessory cart

Accessory Name:
Accessory cart

Application:
To place general accessories for the operating table

Compatible with:

8500/8300
HyBase

6100
HyBase

3000
HyBase

30
UniBase

115-008220-00

Comfort Solution Accessories



- 115-005872-00 Head ring,large, 20×7×7
- 115-005873-00 Head ring,middle, 20×7×5
- 115-005877-00 Head ring,small, 20×7×3
- 115-005874-00 Paediatric head ring,large,14×5×3
- 115-005878-00 Paediatric head ring,middle, 14×7×2.5
- 115-005875-00 Paediatric head ring,small, 10×4×2.5
- 115-005876-00 Infant head ring, 8×3×2



- 115-005885-00 Horseshoe head ring, large, 20×7×7
- 115-005886-00 Horseshoe head ring, middle, 20×7×5
- 115-005887-00 Paediatric horseshoe head ring, L14×5×3
- 115-005888-00 Paediatric horseshoe head ring, S10×4×2.5
- 115-005889-00 Infant head ring, 8×3×2
- 115-005890-00 Horseshoe head ring, small, 20×7×3



- 115-005896-00 Ophthalmic head pad, large, 23×23×4
- 115-005897-00 Ophthalmic head pad, small, 22×19×6
- 115-005898-00 Pillow for supine position, 24×22×7
- 115-038076-00 Ophthalmic pad, paediatric, 15×15×3.5



- 115-005910-00 Lateral positioner, 72×46×15



- 115-005788-00 Pad for intervertebral disc, 63×46×15
- 115-038080-00 Prone position pad, hard inner, 60×49×18cm



- 115-038080-00 Prone position pad, hard inner, 60×49×18cm




- 115-005864-00 Prone position pillow, large, 28×24×14
- 115-005830-00 Prone position pillow, small, 25×20×10
- 115-038074-00 Prone pillow, full support, 23×20×11
- 115-038075-00 Paediatric prone pillow, full 21×18×11



- 115-005892-00 Universal supine pad, 18×18×8



- 115-005900-00 Universal head pad, large, 25×23×6
- 115-005832-00 Universal head pad, small, 25×18×3

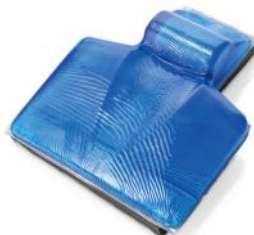

- 115-038077-00 Struma pad, 48×39×9



- 115-005803-00 Short gel body roll, large, 30×15×7
- 115-005807-00 Short gel body roll, middle, 30×10×7
- 115-005804-00 Short gel body roll, small, 30×10×4
- 115-038079-00 Gel body roll, paediatric, 20×10×7cm



- 115-005805-00 Long gel body roll, large, 50×15×7
- 115-005808-00 Long gel body roll, middle, 50×10×7
- 115-005806-00 Long gel body roll, small, 50×10×4

Comfort Solution Accessories

Comfort Solution Accessories



- 115-005796-00 Body strap, 60×7×1
- 115-005906-00 Elbow protection strap, 40×15×1.5
- 115-005822-00 Ankle protection strap, 45×7×1
- 115-005908-00 Wrist protector, 25×5×1



- 115-005790-00 Chest-hipbone pad, 50×20×7



- 115-005794-00 Pad for contraction post, 27×20×1.5



- 115-005814-00 Heel pad, one piece, 18×10×7



- 115-038078-00 Universal hip pad, 50×40×1.5cm



- 115-005826-00 Table pads



- 115-005903-00 Universal arm board pad, short, 40×16×1
- 115-005904-00 Universal arm board pad, long, 60×16×1



- 115-005811-00 Arm/leg pad, long, 50×15×4
- 115-005812-00 Arm/leg pad, short, 40×15×4



- 115-005818-00 Universal square pad-small, 30×20×1.5
- 115-005819-00 Universal square pad-middle, 40×24×1.5
- 115-005820-00 Universal square pad-large, 50×30×1.5
- 115-038078-00 Universal hip pad, 50×40×1.5



- 115-047824-00 Lateral positioner 64×46×25cm

Comfort Solution Accessories