

Radarmed 650+ and 950+



Unit for microwave therapy



Radarmed 650+ and 950+

The Radarmed 650+ and 950+ are units for continuous and pulsed microwave therapy with a frequency of 2450 MHz (12 cm waves). The options for dosage are based on the warmth sensation of the patient. This is considered a big advantage and also makes microwave a safe therapy. Microwave therapy is simple to apply thanks to the unipolar application. The radiator makes no contact with the patient and is simple, safe and hygienic to apply.

- Continuous or pulsed output
- Effortless transition between continuous and pulsed therapy
- The parameters are conveniently arranged on the illuminated display.
- The parameter settings remain visible on the display during treatment.
- Free memory for the storage of treatment protocols (only on Radarmed 950+)

Radarmed 650+ and 950+

1435901 Radarmed 650+

Unit for continuous en pulsed microwave therapy. Dosage and treatment time are digitally displayed on screen. They are easy to set with the push buttons.

1435921 Radarmed 950+

Like the Radarmed 650+, but with LCD-screen. User-friendliness and functionality are essential. The parameters are conveniently arranged on the illuminated display. The parameters are easily and swiftly installed by means of the central control and push buttons. Besides, treatment protocols can be saved.

Standard accessories

- 1435750/51 Operating manual Radarmed 650+/950+(1)
- 1462560 Arm for radiator of the Radarmed
- 1409813 HF-cable 250 Watt, for Radarmed
- 3444298 Mains cable 230V-EUR

ORDERING DATA



650+



950+

(1) = model dependent

Radarmed 650+ and 950+

TECHNICAL SPECIFICATIONS

Frequency	2450 MHz +/- 50 MHz
HF power, continuous	0 - 250 Watt +/- 20% at 50 Ohm
HF power, pulsed	0 - 250 Watt +/- 20% with a fixed peak power of 1500 Watt
Operating modes	continuous and pulsed
Clock	0 - 30 minutes
Mains voltage	220-240 Volt
Frequency	50 / 60 Hz
Current consumption	1100 VA
Safety class	1 type B
Weight	45 kg (including arm)
Dimensions	93x39x49 cm (hwxwd including castors)



Radarmed 650+ and 950+

ACCESSORIES

Optional accessories

- 1409810 Large field radiator for Radarmed
- 1409811 Longitudinal field radiator for Radarmed
- 1409812 Local field radiator Ø 17 cm for Radarmed
- 3441306 Neon checklight
- 1480803 Earth potential-equalization cable 5 m
- 3442022 Protective goggles



1409810



1409811



1409812




3441306



1480803



3442022

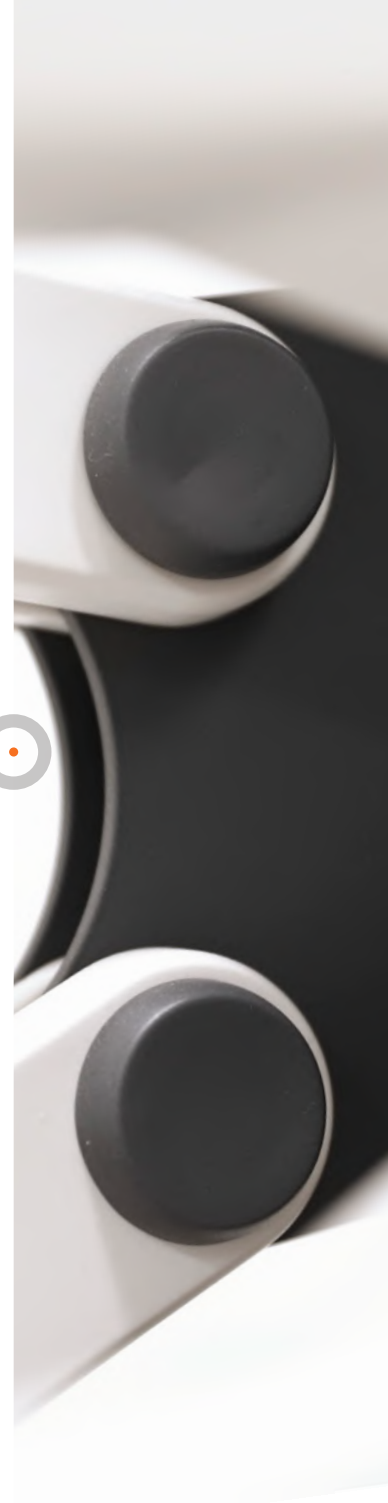

Radarmed 650+ and 950+

Furniture for shortwave and microwave

Special for treatment with shortwave and microwave equipment we offer various furniture. With this furniture there won't occur any heat concentrated areas.

3441401 Treatment chair, clear lacquered wood

3441337 Treatment couch, clear lacquered wood, 190x60x75 cm

3441338 Footrest, clear lacquered wood



ACCESSORIES





Radarmed 650+ and 950+



Unit for microwave therapy

Enraf-Nonius B.V. | Vareseweg 127 | 3047 AT Rotterdam | The Netherlands
www.enraf-nonus.com | info@enraf-nonus.nl | + 31 - (0) 10 20 30 600



PARTNER FOR LIFE